



YOU DRIVE OUR DESIGN

QUICKIE wheelchairs are designed for you and by you. In development, we consider how you interact with your wheelchair to maintain maximum efficiency, performance, fit, and style.

LEADING TECHNOLOGY

We are continually developing new products, improving our flagship products, and offering new options and solutions to enhance your lifestyle.

QUALITY AND DURABILITY

QUICKIE wheelchairs undergo rigorous testing that surpasses industry standards, so you can be confident that your product will hold up and perform throughout its lifetime.


ADJUSTABILITY AND VERSATILITY

Once you choose a QUICKIE, our goal is that it continues to be the right chair for you. We offer maximum adjustability in our designs, interchangeable parts, and the True Fit Growth Program* to help accommodate changes in your condition, environment or lifestyle.

*Available on select QUICKIE wheelchairs

BUILT-4-ME PERSONALISED

Our Built-4-Me Team is dedicated to creating product user-matched beyond what is listed on our product order forms. Experts in the engineering and manufacturing of configured rehab products will listen to your specific requests and build a chair as individual as your fingerprint!

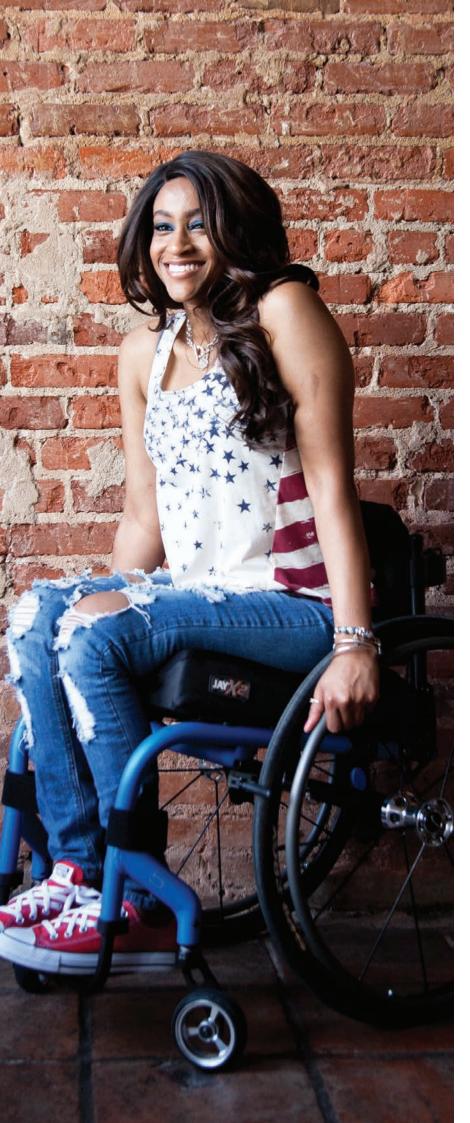



	Live without limits.	
RIGID FRAME		
Krypton R	4	
Nitrum	5	
7 ^R Series	6	
5 ^R	7	
GP Series	8	
QRi	9	
RIGID Comparison (Guide 10	
FOLDING FRA	AME	
Krypton F	13	
Xenon ²	14	
Quickie 2	15	
QXi	16	
M6	17	
FOLDING Comparison	on Guide 18	
TILT IN SPAC	E	
Iris	20	
Iris with Power Tilt	21	
SR45	22	

JAY® Seating	
JAY® Seating Selection	24

TILT IN SPACE Comparison Guide

23

KRYPTON R

CARBON FIBRE FRAME

REVOLUTIONARY

CARBON FRAME

WITH STAGGERING

STRENGTH-TO-WEIGHT RATIO.

RIGID FRAME

Ultra-lightweight from 3.65kgs in transit




CARBON BRAID TECHNOLOGY



This QUICKIE Krypton Range features a patented and fully automated BRAID process that makes it possible to build the perfect frame – strong, light and manufactured with a precise and reliable automated process.

RIGIDITY & RAZOR-SHARP AGILITY

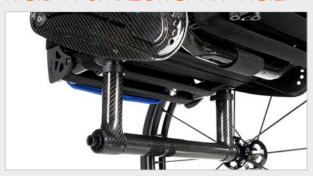

It's not only the exquisite material or the famously seductive styling, but also the athleticism and agility of a unique rigid wheelchair. Infinitely adjustable centre of gravity, fine-tuning of the front and rear seat heights and backrest angle adjustment ensures the Krypton R is designed for something truly special.

NITRUM LIGHTEST IN CLASS

ENGINEERED
TO EFFICIENTLY
CONVERT YOUR
ENERGY INTO MOTION

Available in:

OPEN FRAME DUAL FRAME

Nitrum boasts the largest range of sizes and adjustments in its class, so you can achieve your own perfect fit




LIGHT AND BACK SUPPORT



LED Lights mounted on each caster arm allow you to see and be seen.

One-handed folding back support offers you simple one-handed release and a clean look.

PATENTED, ADJUSTABLE CASTER SYSTEM

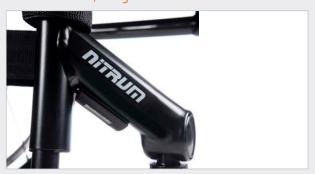

Instead of several welded parts, we use a reinforced and patented caster mechanism in a single hollow- forged component.

Fine-Tuned Adjustment

- Angle
- Camber
- Height

7^R SERIES

RIGID ACTIVE FRAME


STRENGTH IN
MOTION
DESIGNED TO MOVE
WITH YOU!

Available in:

7^R - ADJUSTABLE 7^{RS} - FIXED



FREESTYLE BACKREST

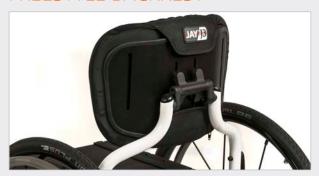

The 7 Series' minimal, ultra-modern design takes backstage so YOU are in the limelight.

Designed to break barriers to independence by providing a full range of upper body motion and lowering the overall weight of your chair.

TOP OF THE LINE TECHNOLOGY



The genius behind the 7R Series' ultra lightweight frame is the 7000 Series aerospace aluminum that has been heat treated using ShapeLoc Technology.

The Result: thin, light tubing with strength & stiffness with a transport weight starting at 4.5kgs*

*Transport weight for the 7RS





HIGHLY
ADJUSTABLE
TO ELEVATE YOUR
MOBILITY

Frame Styles:

CLASSIC

Traditional frame length for maximum stability

PERFORMANCE

Shortened length for maximum responsiveness



AXLE PLATE DESIGN



The 5R's swept-back axle plate design isn't just for looks

It is also incredibly quick to adjust and is available in two lengths to ensure you are getting the adjustments you need without adding extra weight.

XTR SUSPENSION



The XTR Suspension Frame helps you fearlessly conquer curbs, bumps, and tough terrain. Add Mountain wheels, Knobby tires, wide casters, and Frog Legs suspension caster forks for the ultimate off-road wheelchair experience.

GP SERIES

RIGID WHEELCHAIR

A CLASSIC QUICKIE RIGID WHEELCHAIR

Available in:

GP - Standard Style

GPV - Tapered for increased maneuverability

GPSA - Swing Away

for easy transfers



THE ORIGINAL RIGID FRAME



The GP's classic box frame and one piece construction exceed our superior standard for quality and deliver an unmatched rigid feel.

SPORTS OPTION




The tapered GPV front frame, optional Quickie Performance wheels, and up to 12° of wheel camber combine for exceptional maneuverability on and off court.



NO-NONSENSE
UTILITY AND
UNBEATABLE VALUE

COMMUTER FRIENDLY DESIGN

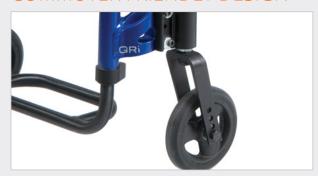


- Ultra lightweight and durable aluminum frame, WC-19 Transit approved up to 120kgs
- Open frame and fluid caster housing design for a stylish look
- Small footprint for maneuverability and easy transport

COMFORTABLE FIT



The QRi comes standard with 0" or 1" of frame inset, a highly adjustable axle plate for perfecting your Center-of-Gravity setting, and multiple front frame angle options.

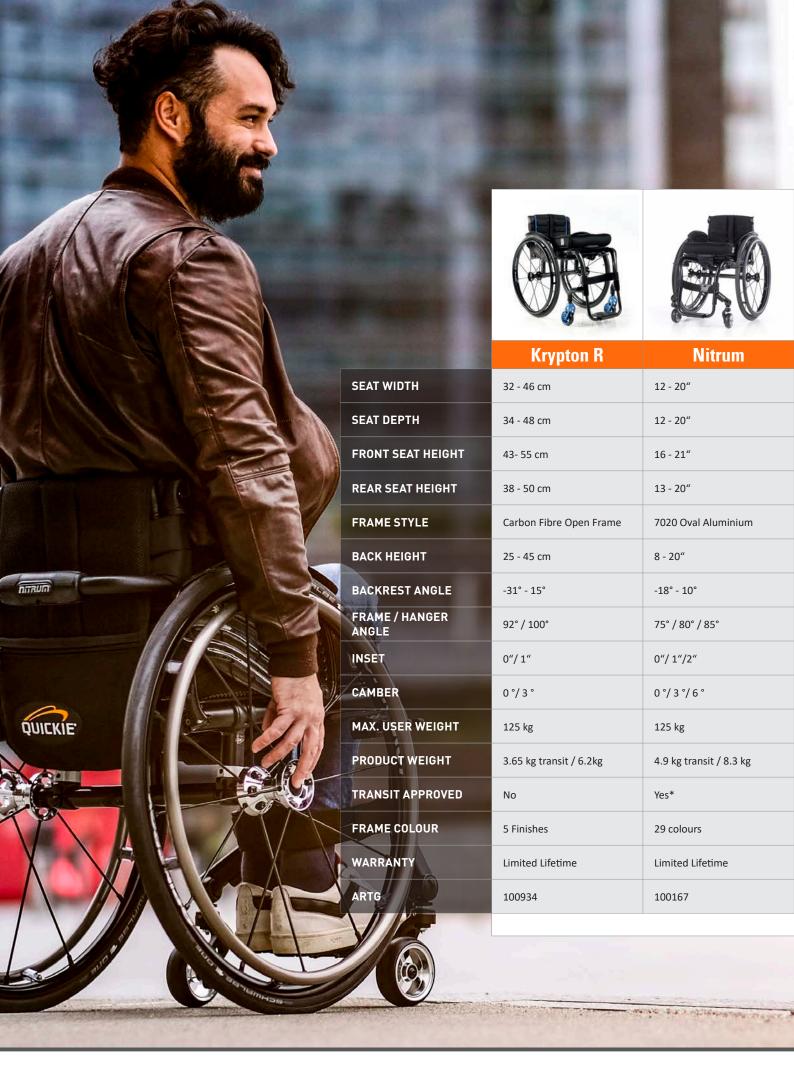

RIGID FRAME SPECIFICATIONS











7 ^R	5 ^R	GP SA	GP/GPV	QRi
12 - 20"	12 - 20"	12 - 20"	12 - 22"	12 - 20"
12 - 20"	12 - 22"	12 - 22"	12 - 22"	12 - 20"
16 - 21"	16 - 21"	16.5 - 23"	15.5 - 23"	15 - 20"
13 - 20"	13 - 20"	16 - 20.5"	14 - 20.5"	14 - 20"
Open Series 7000 Aluminium	Open Series 5000 Aluminium	Traditional	Traditional	Open Series 5000 Aluminium
8 - 20"	10 - 20"	8 - 20"	8 - 20"	8 - 20"
-27° - 8°	-17° - 8°	-4° - 8°	-4° - 8°	-8° - 17°
75° / 80° / 85°	75° / 80° / 85°	60° / 70° / 90°	60° / 70° / 75°/ 80°	75° / 80°/ 85°
0"/ 1"/2"	0"/ 1"/2"	N/A	N/A on GP, 1" V taper on GPV	0, 1" V inset
0°/3°/6°	0°, 3°, 6°	(Adjustable) Set at 3°	(Adjustable) Set at 3°	0°/3°
120 kg	120 kg	113 kg	113 kg	120 kg
5.1 kg transit / 8.5 kg	6.4 kg transit / 9.8 kg	8.6 kg transit / 12 kg	7.6 kg transit / 11 kg GP	6.6 kg transit / 10 kg
Yes*	Yes*	ISO Only - thru frame	N/A on GP, ISO Only - thru frame on GPV	Yes*
29 colours	28 colours	28 colours	28 colours	24 colours
Limited Lifetime	Limited Lifetime	Limited Lifetime	Limited Lifetime	Limited Lifetime
100167	100167	100167	100167	100167

 $* Conforms \ to \ ANDI/RESNA \ WC-19 \ Standards. \ Please \ refer \ to \ the \ Quickie \ \& \ Zippie \ WC-19 \ Transit \ Securement \ Supplement$





KRYPTON F CARBON FIBRE FRAME

REVOLUTIONARY
CARBON FRAME
WITH STAGGERING
STRENGTH-TOWEIGHT RATIO

FOLDING FRAME

Ingeniously designed crossbrace with an open frame feel!



CARBON BRAID TECHNOLOGY



This QUICKIE Krypton Range features a patented and fully automated BRAID process that makes it possible to build the perfect frame – strong, light and manufactured with a precise and reliable automated process.

LIGHTEST FOLDING WHEELCHAIR



From just 5.8kgs in transit, the QUICKIE Krypton F is our lightest folding wheelchair yet! From every angle, the Krypton F is conceived to portray beauty. The heart of the Krypton F is its ingeniously designed cross-brace that fits so neatly under the seat it's barely noticeable giving a very minimalistic, open-frame feel.

GUICKIE 2 FOLDING WHEELCHAIR

ONE OF THE LIGHTEST FOLDING

WHEELCHAIRS IN

THE WORLD!

Available in:

FF - Open fixed front frame


DUAL - Reinforced

fixed front frame

SA - Swing Away for easy transfers



EASY TO FOLD MECHANISM

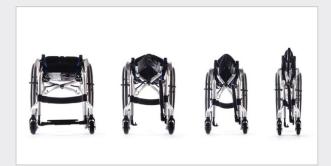

The aluminium latches on the cross-brace will also keep your wheels in parallel when folding and makes sure the chair stays folded when lifting.

ADJUSTABLE TO SUIT YOU



You want to change your chair driving characteristic from passive to more active or the opposite way? Just change the centre of gravity setting by simply repositioning the axle stem. And don't forget to adjust the back angle for your ergonomic posture.



VERSATILE AND MODULAR.
THE INDUSTRY'S

FAVORITE FOR OVER

25 YEARS!

Available in:

Q2 - Standard Modular frame

Q2 LITE - *Ultra lightweight, modular frame*

Q2 HP - High Performance



MODULAR FRAME DESIGN



Featuring a modular design, the QUICKIE 2 offers the broadest range of possible configurations in the industry: over 100,000 different complex rehab options and counting. Designed to meet all of your unique needs without compromising fit or performance, only one wheelchair can live up to the name QUICKIE.

QUICKIE 2 LITE



Love the Quickie 2, but looking for something with an edge? The Quickie 2's minimalist counterpart, the Quickie 2 Lite, features an extremely lightweight, modular frame with a rounded rear for unique style and either fixed or swing-away front hangers.

FOLDING WHEELCHAIR

LOW
MAINTENANCE
& SIMPLICITY



DURABLE FRAME DESIGN




- Ultra lightweight and durable folding frame with a 120kgs weight capacity
- One-piece side frames reduce moving parts
- Three arm cross brace sits inside frame for rigid ride

SIMPLE ADJUSTMENT & SERVICING



Most adjustments can be made in less than one minute, using only one tool - Center-of-Gravity (axle plate), Seat-to-floor height, Camber, Axle spacing and Back height.



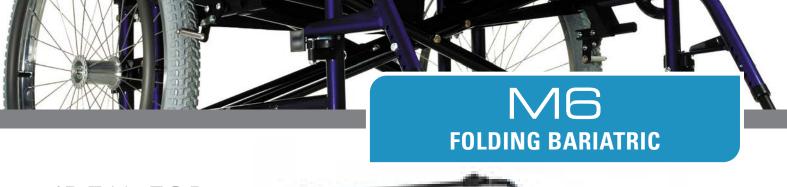



MODULAR & HIGH-STRENGTH



The M6's modular frame, features high strength steel for increased durability. Its interchangeable, heavy duty components accommodate changing conditions with simple parts replacements.

LOCKING CROSS BRACE & STABILSER BAR



The M6 utilises a double locking cross brace and stabiliser bar to minimize frame flex and increase durability. The frame will hold an active user upto 295 kg.





FOLDING FRAME SPECIFICATIONS











Xenon² Dual®	Quickie 2 Lite	Quickie 2	QXi	M6
32 - 50 cm	14 - 20"	12 - 24"	14 - 20"	22 - 30"
34 - 50 cm	14 - 20"	12 - 22"	14 - 20"	18 - 22"
43 - 55 cm	14.5 - 21.5"	13.5 to 22"	12.5 - 19.5"	17 - 20"
37 - 50 cm	13 - 19.5"	13 to 22"	12 - 19"	17 - 20"
Reinforced with fixed front	Traditional - Fixed or SA front	Traditional- Fixed or SA front	Traditional	Heavy Duty
25 - 47.5 cm	9 - 21"	9 - 21"	13 - 20"	17 - 20"
75° to 103° in 4° increments	-14° - +5° (optional)	-15° - +3° (Optional)	-15° - 3°	N/A
92° / 100°	70°/ 80°	60°/ 70°/ 80°/ 90°	60°/ 70° / 80°	N/A
0 - 1 cm	0 / 1" each side	N/A	N/A	N/A
0° / 2° / 4°	(Adjustable) Set at 3°	(Adjustable) Set at 3°	(Adjustable) Set at 3°	0°
140 kg	120 kg	136 kg Std / 158 kg HD	136 kg	295 kg
7.2 kg transit / 9.5 kg	7.3 kg transit / 10.6 kg	8.2 kg transit / 11.5 kg	9.9 kg transit / 13.25 kg	24 kg
Yes	Yes	Yes	Yes	No
32 colours	24 colours	24 colours	24 colours	24 colours
Limited Lifetime	Limited Lifetime	Limited Lifetime	Limited Lifetime	Limited Lifetime
100934	100167	100167	100167	100167

IRIS TILT IN SPACE

INTELLIGENT ROTATION IN SPACE TECHNOLOGY

One of the best and lightest tilt-in-space wheelchairs on the market



INTELLIGENT ROTATION



The Quickie IRIS features up to 55° of Intelligent Rotation in Space Technology to achieve superior positioning and is designed to facilitate feeding and respiratory function, reduce pressures and improve visual alignment.

MONO BACKREST SYSTEM




The MONO Backrest System's single back post offers limitless seating and positioning capabilities. By mounting the backrest to a single central post, the backrest and chair widths may be selected independently to accommodate users with varying torso and hip widths.

IRIS WITH POWER TILT

HAVE MORE
CONTROL OVER THE

POSITION OF YOUR

BODY IN SPACE

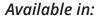

IRIS

IRIS Heavy Duty

- improves comfort by assisting with sitting tolerance
- provides more independence to participate in daily activities



ATTENDANT CONTROL



Allows the carer to give control of the tilt function to the user. The battery gauge is displayed on the top of the control so you never have to guess when it's time to charge. Ability to conveniently hang the attendant control anywhere on the chair.

TOGGLE SWITCH



Can be mounted vertically or horizontally depending on needs. Charging port is conveniently located on the toggle switch for easy access.


SR45 TILT IN SPACE

PRE-CONFIGURED

FOR EASE OF ORDERING



45° OF ROTATION-IN-SPACE

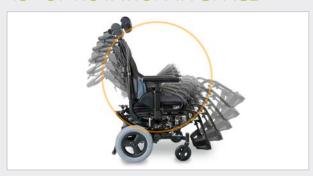


The Quickie SR45's seat frame rotates around the client's centre of gravity providing the shortest possible wheelbase without compromising stability. When the centre of gravity and centre of rotation are properly aligned, tilting is smooth, easy and controlled.

STANDARD SMART PACKAGE



SR45 comes standard with - Non folding back canes, Height adjustable back with adjustable stroller handles, Foot tilt, Transit tie downs, Dual post height adjustable flip back with full length arm pads, 70` Swing in/ out hangers, solid seat pan and rear anti tips and composite foot plates.



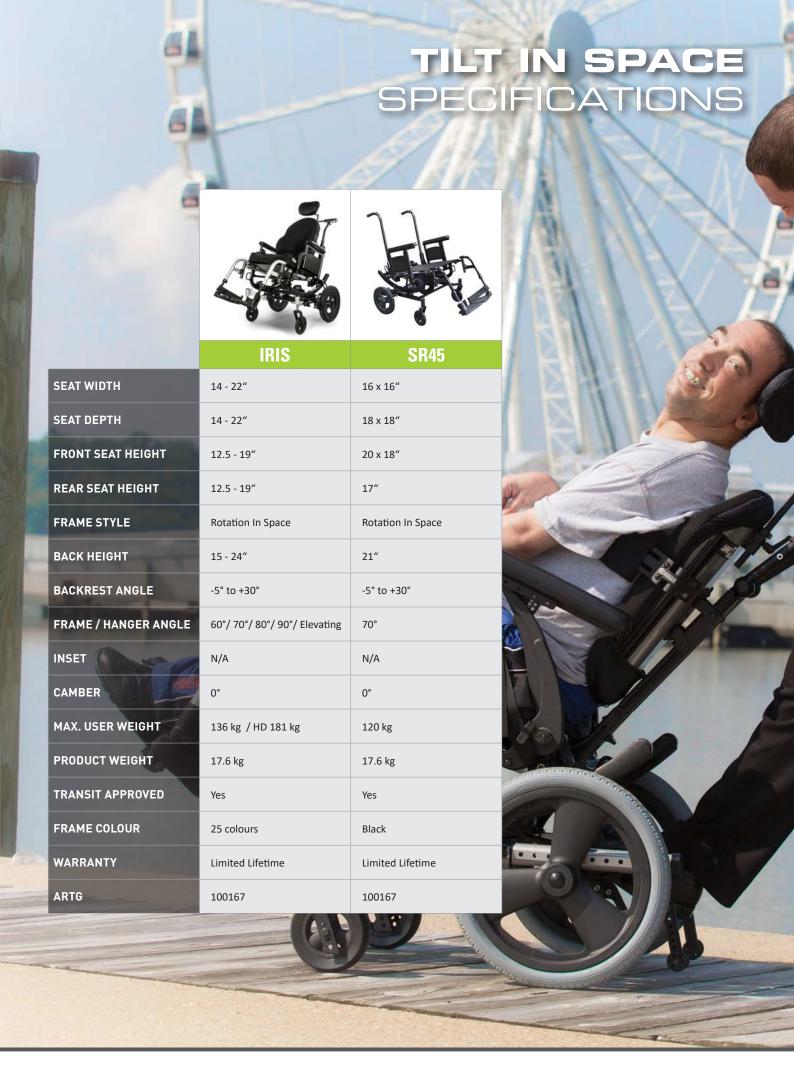





SEATING

JAY® BACKS

JAY® CUSHIONS



Basic

Basic PRO



 $G0^{\scriptscriptstyle \mathsf{TM}}$





Soft Combi™ P



Basic





Easy™



Ion™







J3™ Plus **Focus Point**













J3™ Posterior

J3™ Posterior & Lateral

J3™ PD, PDL & PDC

To learn about personalising JAY seating products, visit the JAY Your Way product page on www.SunriseMedical.com.au



Sunrise Medical Australia Pty Ltd 11 Daniel Street, Wetherill Park, NSW 2164 Phone: (02) 9678 6600 Email: enquiries@sunrisemedical.com.au Website: www.sunrisemedical.com.au